


August 5, 2015

As required for compliance to the latest ISO Standard (ISO 81060-1), components of blood pressure measurement systems are no longer allowed to be supplied with luer-lock or luer-slip compatible connectors. As a result, the following product changes have been implemented on blood pressure cuffs and aneroid gauges.

SCREW Configurations:**FlexiPort -1SC & -2SC configurations****Trimline 1 & 2-tube screw configurations (including Neonatal cuffs) Private Labeled 1 and 2-tube screw configurations:**

1. The female screw connector used on the above referenced cuffs traditionally had an internal lumen profile that would allow insertion and capture of a male luer slip connector. We have eliminated that feature in all Screw type connectors to be compliant to ISO 81060. This has not affected the ability to mate with device side (male) screw connectors.

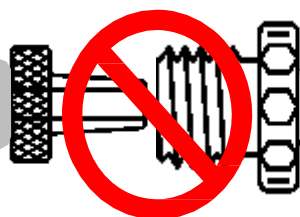

2. Trimline disposable and reusable cuffs with female luer slip connectors have been discontinued.
3. Neonatal Cuffs with male luer connectors remain available in the US and other non-CE countries.

What do Customers need to do:

If a customer was using the cuff side screw connector as a luer slip, they will need to change the connectors on their Device hoses to the mated-screw connector which is part # **5082-161** (10 per bag). If the facility/hospital has many devices that need to be changed, they can contact their Field Sales Representative to discuss other options, for example, converting to the FlexiPort connector which is compliant to the ISO 81060 standard and Welch Allyn's conversion team can assist them with the change of their connectors.

TriPurpose (TP) Connector:

The Female Luer Slip and Female Luer Lock capabilities of the TriPurpose connector are being removed in two phases. The first phase of the change resulted in a reduction to the inside diameter of the TP connector which now prevents insertion of male luer-slip and luer-lock connectors. Once tooling changes are complete, the outside male threads that comprise the “lock” feature of the female luer lock will also be removed. The TriPurpose connector is now a “UniPurpose” connector, but the TP code in all part numbers will remain unchanged.



Phase 1 (Completed): Cuff side TP Connector looks unchanged but female luer profile has been removed.



Phase 2: Once tooling changes are complete, Male Luer Lock thread will be removed leaving only the TP Male screw that mates with the device side female TP connector

What do Customers need to do:

Customers that were using the tripurpose cuff side connector as a luer-lock or luer slip can change their device side connector to the female Tri-Purpose connector which is part # **5082-100** (10 per bag). The 5082-100 is compatible with the redesign of the cuff side tri-purpose connector.

DS65 & DS66 (Trigger Style) Aneroid Gauge:

Welch Allyn has replaced the Male Luer Lock connector with the Female Tri-Purpose connector on all DS65 and DS66 Trigger Style Aneroid Gauges



Before
Luer-Lock Male



After
TP Female

What do Customers need to do:

Customers that have devices with the luer-lock fitting can change their device side connector to the female Tri-Purpose connector which is part # **5082-100** (10 per bag).

For those interested in the specific clauses in the new ISO standard, these are listed below:

Subclause 7.2.6 Tubing connectors

This requirement addresses the known hazards associated with cross-connections. There is a demonstrated unacceptable risk associated with the use of Luer connectors with patients, including neonatal patients, in the cuffs and tubing. Three-way adaptors are readily available in the clinical environment, so that a reverse Luer connector affords very little protection from this risk.

7.2.6 Tubing connectors

Tubing connectors, if provided, shall incorporate a means of preventing accidental disconnection. Tubing connectors shall not be equipped with a connector that connects with a connector complying with ISO 594-1 or ISO 594-2.

Normative references

ISO 594-1, *Conical fittings with a 6 % (Luer) taper for syringes, needles and certain other medical equipment — Part 1: General requirements*

ISO 594-2, *Conical fittings with 6 % (Luer) taper for syringes, needles and certain other medical equipment —*

Regards,

David Tufenkjian
Global Category Specialist
Core Blood Pressure
315 685-4115

Attachment: List of affected BP cuff part numbers.

Material #	Material Description	Cuff Style	NOTES ON CHANGE
1360	Connector, 5/32" Hose Barb to FEMALE Threaded, 10/PK	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
1376	Connector, 5/32" Hose Barb to MALE Threaded, 10/PK	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9001	ADAPTER,DBL SC to SGL 1/8in BARB	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9002	ADAPTER,DBL SC to DBL 1/8in BARB	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9003	ADAPTER,DBL SC to SGL ML	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9072	ADAPTER,SGL SC to DBL 1/8in BARB	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9073	ADAPTER,SGL SC to SGL ML	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9074	ADAPTER,SGL SC to DBL ML	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9097	ADAPTER,DBL MF to DBL MSC	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9206	ADAPTER,DBL SC to DBL ML	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9218	ADAPTER,DBL MSC to DBL ML	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9229	ADAPTER,DBL SC to DBL FL	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9230	ADAPTER,SGL HP to DBL MSC	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9231	ADAPTER,SGL LL to DBL MSC	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9236	ADAPTER,DBL MQ to DBL MSC	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9240	ADAPTER,DBL SC to DBL NK	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9249	ADAPTER,SGL 5/32in BARB to DBL MSC	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9268	ADAPTER,SGL SC to DBL NK	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9269	ADAPTER,SGL SC to DBL MSC	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9270	Adapter, Double 5/32" Hose Barb & Bare Tube to Double Male Threaded Connector, 1/EA	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9272	ADAPTER,SGL EXT TUBING to DBL MSC and BV	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9274	ADAPTER,DBL FL and BV to DBL MSC	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161

Material #	Material Description	Cuff Style	NOTES ON CHANGE
9281	ADAPTER,DBL SC to DBL FL	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9282	ADAPTER,DBL SC to SGL 5/32in BARB	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9286	ADAPTER,SGL SC to SGL FHP METAL	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9288	ADAPTER,DBL SC to SGL FHP METAL	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9290	ADAPTER,SGL HP to SGL MSC	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9291	ADAPTER,SGL MQ to DBL MSC	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9292	ADAPTER,DBL FL 5/32in BARB to DBL MSC	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9293	ADAPTER,SGL ML to DBL MSC	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9294	ADAPTER,DBL BV and EXT TUBING to DBL MSC	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9298	ADAPTER,DBL ML to DBL MSC	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9300	ADAPTER,SGL NK to SGL MSC	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9302	ADAPTER,SGL NK to DBL MSC	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9307	ADAPTER,DBL SC to SGL MSC	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9308	ADAPTER,SGL SC to DBL MF	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9314	ADAPTER,DBL SC to SGL FHP	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9319	ADAPTER,DBL SC to SGL MQ	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9329	ADAPTER,SGL SC to SGL FHP	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9334	ADAPTER,DBL FL to DBL MSC	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9336	ADAPTER,DBL SC to DBL MF	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9351	ADAPTER,SGL SC to DBL FHP	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9352	ADAPTER,DBL SC to SGL FL	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9354	ADAPTER,DBL FL and 10in BT to DBL MSC	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161

[illegible]

[illegible]

Material #	Material Description	Cuff Style	NOTES ON CHANGE
39218	CUFF TL DISP,SOFT,THIGH,2TUBE,SCC	Trimline Soft	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
39229	CUFF TL DISP,SOFT,LG AD LONG,2TB,BVSCC	Trimline Soft	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
39235	CUFF TL DISP,SOFT,SM CHILD,2TUBE,BVSCC	Trimline Soft	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
39257	CUFF TL DISP,SOFT,CHILD LONG,2TB,BVSCC	Trimline Soft	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
66018	Connector, 1/8" Hose Barb to MALE Threaded/ 10/PK	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
66024	Connector, 3/16" Hose Barb to MALE Threaded, 10/PK	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
66027	ADAPTER,SGL LL to SGL MSC	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
66033	ADAPTER,SGL MLL to SGL MSC	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
66038	ADAPTER,SGL BT to SGL MSC	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
66041	ADAPTER,SGL HP to SGL MSC	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
66044	ADAPTER,SGL LL to DBL MSC	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
66047	ADAPTER,SGL MLL to DBL MSC	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
66048	ADAPTER,SGL SC to DBL MSC	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
66050	ADAPTER,SGL BT to DBL MSC	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
66053	ADAPTER,SGL HP to DBL MSC	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
66158	ADAPTER,SGL MQ to SGL MSC	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
99702	CUFF TL DISP,VINYL,NEO2,2TUBE,SC	Neonatal Vinyl TL	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
99715	CUFF TL DISP,VINYL,NEO3,1TUBE,SC	Neonatal Vinyl TL	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
99724	CUFF TL DISP,VINYL,NEO1,1TUBE,SC	Neonatal Vinyl TL	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
99748	CUFF TL DISP,VINYL,NEO5,1TUBE,SC	Neonatal Vinyl TL	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
99755	CUFF TL DISP,VINYL,NEO3,2TUBE,SC	Neonatal Vinyl TL	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
99761	CUFF TL DISP,VINYL,NEO4,2TUBE,SC	Neonatal Vinyl TL	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161

[illegible]

Material #	Material Description	Cuff Style	NOTES ON CHANGE
99878	CUFF TL DISP,VINYL,SM ADULT,2TUBE,SCC	Trimline VINYL	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
99879	CUFF TL DISP,VINYL,ADULT,2TUBE,SCC	Trimline VINYL	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
99880	CUFF TL DISP,VINYL,ADULT LONG,2TB,SCC	Trimline VINYL	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
99881	CUFF TL DISP,VINYL,LG ADULT,2TUBE,SCC	Trimline VINYL	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
99882	CUFF TL DISP,VINYL,ADULT,2TUBE,BVSCC	Trimline VINYL	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
99883	CUFF TL DISP,VINYL,LG AD LONG,2TB,SCC	Trimline VINYL	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
99884	CUFF TL DISP,VINYL,THIGH,2TUBE,SCC	Trimline VINYL	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
99885	CUFF TL DISP,VINYL,LG ADULT,2TUBE,BVSCC	Trimline VINYL	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
99886	CUFF TL DISP,VINYL,ADULT LONG,2TB,BVSCC	Trimline VINYL	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
1360N	Connector, 1/8" Hose Barb to FEMALE Neonatal Threaded, 10/PK	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
1900CK	CUFF TL REUSE,1 PC,THIGH,2TUBE,SC	Trimline One-Piece Reuse	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
1900SCK	CUFF TL REUSE,1 PC,THIGH,1TUBE,SC	Trimline One-Piece Reuse	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
1902CK	CUFF TL REUSE,1 PC,LG ADULT,2TUBE,SC	Trimline One-Piece Reuse	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
1902SCK	CUFF TL REUSE,1 PC,LG ADULT,1TUBE,SC	Trimline One-Piece Reuse	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
1902XLCK	CUFF TL REUSE,1 PC,LG AD LONG,2TB,SC	Trimline One-Piece Reuse	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
1902XLCK	CUFF TL REUSE,1 PC,LG AD LONG,1TB,SC	Trimline One-Piece Reuse	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
1903CK	CUFF TL REUSE,1 PC,ADULT,2TUBE,SC	Trimline One-Piece Reuse	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
1903SCK	CUFF TL REUSE,1 PC,ADULT,1TUBE,SC	Trimline One-Piece Reuse	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
1903XLCK	CUFF TL REUSE,1 PC,ADULT LONG,2TB,SC	Trimline One-Piece Reuse	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
1903XLCK	CUFF TL REUSE,1 PC,ADULT LONG,1TB,SC	Trimline One-Piece Reuse	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
1904CK	CUFF TL REUSE,1 PC,SM ADULT,2TUBE,SC	Trimline One-Piece Reuse	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
1904SCK	CUFF TL REUSE,1 PC,SM ADULT,1TUBE,SC	Trimline One-Piece Reuse	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161

Material #	Material Description	Cuff Style	NOTES ON CHANGE
1905CK	CUFF TL REUSE,1 PC,CHILD,2TUBE,SC	Trimline One-Piece Reuse	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
1905SCK	CUFF TL REUSE,1 PC,CHILD,1TUBE,SC	Trimline One-Piece Reuse	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
1906CK	CUFF TL REUSE,1 PC,INFANT,2TUBE,SC	Trimline One-Piece Reuse	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
1906SCK	CUFF TL REUSE,1 PC,INFANT,1TUBE,SC	Trimline One-Piece Reuse	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
1-SC	BP PORT FITTING,1-TUBE,SCREW	FlexiPort Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
2-SC	BP PORT FITTING,2-TUBE,SCREW	FlexiPort Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
2-SCC	BP PORT FITTING,2-TB,SC & CLAMP	FlexiPort Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30503009SA	CUFF, SOFT, SM INFANT, 1-TUBE, SC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30503-010A	CUFF, SOFT, INFANT, 2-TUBE, SC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30503010SA	CUFF, SOFT, INFANT, 1-TUBE, SC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30503-011A	CUFF, SOFT, CHILD, 2-TUBE, SC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30503011SA	CUFF, SOFT, CHILD, 1-TUBE, SC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30503-012A	CUFF, SOFT, SM AD, 2-TUBE, SC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30503012SA	CUFF, SOFT, SM AD, 1-TUBE, SC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30503-013A	CUFF, SOFT, ADULT, 2-TUBE, SC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30503013AA	CUFF, SOFT, AD LONG 2-TUBE, SC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30503013SA	CUFF, SOFT, ADULT, 1-TUBE, SC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30503-014A	CUFF, SOFT, LG AD, 2-TUBE, SC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30503014AA	CUFF SOFT LG AD LONG 1-TUBE SC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30503014LA	CUFF SOFT LG AD LONG 2-TUBE SC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30503014SA	CUFF, SOFT LG AD, 1-TUBE, SC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30503-015A	CUFF, SOFT, THIGH, 2-TUBE, SC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161

Material #	Material Description	Cuff Style	NOTES ON CHANGE
30503015SA	CUFF, SOFT, THIGH, 1-TUBE, SC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30503-019A	CUFF, SOFT SM CHILD 2-TUBE, SC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30503019SA	CUFF, SOFT SM CHILD 1-TUBE, SC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30503-710A	CUFF, SOFT INFANT, 2-TUBE, SCC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30503-711A	CUFF, SOFT, CHILD, 2-TUBE, SCC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30503-712A	CUFF, SOFT, SM AD, 2-TUBE, SCC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30503-713A	CUFF, SOFT, ADULT, 2-TUBE, SCC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30503-714A	CUFF, SOFT, LG AD, 2-TUBE, SCC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30503LSA	CUFF, SOFT, AD LONG 1-TUBE, SC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30506010A	CUFF, VINYL INFANT, 2-TUBE, SC	FlexiPort VINYL	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30506-011A	CUFF, VINYL CHILD, 2-TUBE, SC	FlexiPort VINYL	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30506011SA	CUFF, VINYL CHILD, 1-TUBE, SC	FlexiPort VINYL	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30506-012A	CUFF, VINYL, SM AD, 2-TUBE, SC	FlexiPort VINYL	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30506012SA	CUFF, VINYL, SM AD, 1-TUBE, SC	FlexiPort VINYL	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30506-013A	CUFF, VINYL ADULT, 2-TUBE, SC	FlexiPort VINYL	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30506013SA	CUFF, VINYL ADULT, 1-TUBE, SC	FlexiPort VINYL	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30506014A	CUFF, VINYL LG AD, 2-TUBE, SC	FlexiPort VINYL	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30506014SA	CUFF, VINYL LG AD, 1-TUBE, SC	FlexiPort VINYL	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30506-015A	CUFF, VINYL THIGH, 2-TUBE, SC	FlexiPort VINYL	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30506714LA	CUFF VINL LG AD LONG 2TUBE SCC	FlexiPort VINYL	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
30506A013L	CUFF, VINYL AD LONG 2-TUBE, SC	FlexiPort VINYL	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
3200CK	CUFF TL REUSE,2 PC,THIGH,2TUBE,SC	Trimline Two-Piece Reuse	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161

[illegible]

Material #	Material Description	Cuff Style	NOTES ON CHANGE
3603SCK	CUFF TL REUSE,1 PC,ADULT,1TUBE,SC	Trimline One-Piece Reuse	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
3603XLCK	CUFF TL REUSE,1 PC,ADULT LONG,2TB,SC	Trimline One-Piece Reuse	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
3603XLCK	CUFF TL REUSE,1 PC,ADULT LONG,1TB,SC	Trimline One-Piece Reuse	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
3604CK	CUFF TL REUSE,1 PC,SM ADULT,2TUBE,SC	Trimline One-Piece Reuse	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
3604SCK	CUFF TL REUSE,1 PC,SM ADULT,1TUBE,SC	Trimline One-Piece Reuse	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
3605CK	CUFF TL REUSE,1 PC,CHILD,2TUBE,SC	Trimline One-Piece Reuse	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
3605SCK	CUFF TL REUSE,1 PC,CHILD,1TUBE,SC	Trimline One-Piece Reuse	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
3606CK	CUFF TL REUSE,1 PC,INFANT,2TUBE,SC	Trimline One-Piece Reuse	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
3606SCK	CUFF TL REUSE,1 PC,INFANT,1TUBE,SC	Trimline One-Piece Reuse	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
3655CK	CUFF TL REUSE,HEALTHCARE SET,SC	Trimline One-Piece Reuse	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
3657CK	CUFF TL REUSE,PEDIATRIC SET,SC	Trimline One-Piece Reuse	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
3659CK	CUFF TL REUSE,PHYSICIAN SET,SC	Trimline One-Piece Reuse	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
5200-01	CUFF, ADULT, W/6 IN TUBE	Two-Piece	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
5200-02	CUFF, LARGE ADULT, W/6 IN TUBE	Two-Piece	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
5200-03	CUFF,CHILD PRINT,W/6 IN TUBE	Two-Piece	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
5200-04	BLADDER,ADULT W/6 IN TUBE	Two-Piece	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
5200-05	BLADDER LG.AD. W/6IN TUBE	Two-Piece	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
5200-06	BLADDER,CHILD,W/6 IN TUBE	Two-Piece	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
5200-10	CUFF & BAG, THIGH, W/6 IN TUBE	Two-Piece	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
5200-11	BLADDER,THIGH W/6 INCH TUBE	Two-Piece	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
9274B	ADAPTER,DBL BARB and BV to DBL MSC	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
A39010	AMERINET CUFF,SOFT,INFANT,2TB,SC	Trimline Soft	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161

[illegible]

Material #	Material Description	Cuff Style	NOTES ON CHANGE
M39023	MEDICHOICE,SOFT,ADULT,2TUBE,SC	Trimline Soft	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
M39053	MEDICHOICE,SOFT,ADULT LONG,2TUBE,SC	Trimline Soft	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
M39071	MEDICHOICE,SOFT,LG ADULT LONG,2TUBE,SC	Trimline Soft	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
M39073	MEDICHOICE,SOFT,ADULT,1TUBE,SC	Trimline Soft	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
M39074	MEDICHOICE,SOFT,LG ADULT,1TUBE,SC	Trimline Soft	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
M39075	MEDICHOICE,SOFT,LG ADULT LONG,1TUBE,SC	Trimline Soft	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
M39091	MEDICHOICE,SOFT,CHILD,1TUBE,SC	Trimline Soft	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
M39092	MEDICHOICE,SOFT,SM ADULT,1TUBE,SC	Trimline Soft	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
M39093	MEDICHOICE,SOFT,ADULT LONG,1TUBE,SC	Trimline Soft	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
M9097	MEDICHOICE ADAPTER,DBL MF to DBL MSC	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
M9230	MEDICHOICE ADAPTER,SGL HP to DBL MSC	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
M9231	MEDICHOICE ADAPTER,SGL HP to DBL MSC	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
M9236	MEDICHOICE ADAPTER,DBL MQ to DBL MSC	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
M9269	MEDICHOICE ADAPTER,SGL SC to DBL MSC	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
M9314	MEDICHOICE ADAPTER,DBL SC to SGL FHP	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
M9329	MEDICHOICE ADAPTER,SGL SC to SGL FHP	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
R39012	DISP BP CUFF, SMALL ADULT, 2TUBE, SC	Trimline SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
R39023	DISP BP CUFF, ADULT, 2TUBE, SC	Trimline SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
R39053	DISP BP CUFF, ADULT LONG, 2TUBE, SC	Trimline SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
REUSE-06-1SC	CUFF, REUS, SM INFANT, 1-TUBE, SC	FlexiPort Reuse	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
REUSE-07-1SC	CUFF, REUS, INFANT, 1-TUBE, SC	FlexiPort Reuse	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
REUSE-07-2SC	CUFF, REUS, INFANT, 2-TUBE, SC	FlexiPort Reuse	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161

[illegible]

Material #	Material Description	Cuff Style	NOTES ON CHANGE
SOFT-09-2SC	CUFF, SOFT, CHILD, 2-TUBE, SC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
SOFT-09-2SCC	CUFF, SOFT, CHILD, 2-TUBE, SCC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
SOFT-10-1SC	CUFF, SOFT, SM AD, 1-TUBE, SC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
SOFT-10-2SC	CUFF, SOFT, SM AD, 2-TUBE, SC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
SOFT-10-2SCC	CUFF SOFT SM ADULT 2-TUBE, SCC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
SOFT-11-1SC	CUFF, SOFT, ADULT, 1-TUBE, SC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
SOFT-11-2SC	CUFF, SOFT, ADULT, 2-TUBE, SC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
SOFT-11-2SCC	CUFF, SOFT, ADULT, 2-TUBE, SCC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
SOFT-11L-1SC	CUFF, SOFT, AD LONG 1-TUBE, SC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
SOFT-11L-2SC	CUFF, SOFT, AD LONG 2-TUBE, SC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
SOFT-12-1SC	CUFF, SOFT, LG AD, 1-TUBE, SC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
SOFT-12-2SC	CUFF, SOFT, LG AD, 2-TUBE, SC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
SOFT-12-2SCC	CUFF, SOFT, LG AD, 2-TUBE, SCC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
SOFT-12L-1SC	CUFF SOFT LG AD LONG 1-TUBE SC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
SOFT-12L-2SC	CUFF SOFT LG AD LONG 2-TUBE SC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
SOFT-13-1SC	CUFF, SOFT, THIGH, 1-TUBE, SC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
SOFT-13-2SC	CUFF, SOFT, THIGH, 2-TUBE, SC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
SOFT-13-2SCC	CUFF, SOFT, THIGH, 2-TUBE, SCC	FlexiPort SOFT	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
TL9271	ADAPTER,SGL FL to DBL MSC and BV	Accessory	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
V39011	NOVAPLUS CUFF,SOFT,CHILD,2TUBE,SC	Trimline Soft	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
V39012	NOVAPLUS CUFF,SOFT,SM ADULT,2TUBE,SC	Trimline Soft	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
V39014	NOVAPLUS CUFF,SOFT,LG ADULT,2TUBE,SC	Trimline Soft	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161

[illegible]

Material #	Material Description	Cuff Style	NOTES ON CHANGE
VINYL-12L-2SC	CUFF VINYL LG AD LONG 2-TUBE SC	FlexiPort VINYL	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
VINYL-13-1SC	CUFF, VINYL THIGH, 1-TUBE, SC	FlexiPort VINYL	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
VINYL-13-2SC	CUFF, VINYL THIGH, 2-TUBE, SC	FlexiPort VINYL	Screw connector can no longer function as Luer-Slip, will only function as a threaded screw connector that mates with the device side screw, part # 5082-161
1706	BLADDER ONLY,TL BALANCED,INFANT,2TUBE	Trimline Two-Piece Reuse	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
2700	BLADDER ONLY,TL BALANCED,THIGH,2TUBE	Trimline Two-Piece Reuse	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
2702	BLADDER ONLY,TL BALANCED,LG ADULT,2TUBE	Trimline Two-Piece Reuse	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
2703	BLADDER ONLY,TL BALANCED,ADULT,2TUBE	Trimline Two-Piece Reuse	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
2704	BLADDER ONLY,TL BALANCED,SM ADULT,2TUBE	Trimline Two-Piece Reuse	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
2705	BLADDER ONLY,TL BALANCED,CHILD,2TUBE	Trimline Two-Piece Reuse	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
1706S	BLADDER ONLY,TL BALANCED,INFANT,1TUBE	Trimline Two-Piece Reuse	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
1-TPL	BP PORT FITTING,1TUBE,TRI-PURP	FlexiPort Accessory	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
1-TPLE	BP PORT FITTING,1-TB,LG EURO	FlexiPort Accessory	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
1-TPS	BP PORT FITTING,1TUBE,TRI-PURP	FlexiPort Accessory	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
1-TPSE	BP PORT FITTING,1-TB,SM EURO	FlexiPort Accessory	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
2700S	BLADDER ONLY,TL BALANCED,THIGH,1TUBE	Trimline Two-Piece Reuse	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
2702S	BLADDER ONLY,TL BALANCED,LG ADULT,1TUBE	Trimline Two-Piece Reuse	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
2703S	BLADDER ONLY,TL BALANCED,ADULT,1TUBE	Trimline Two-Piece Reuse	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
2704S	BLADDER ONLY,TL BALANCED,SM ADULT,1TUBE	Trimline Two-Piece Reuse	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
2705S	BLADDER ONLY,TL BALANCED,CHILD,1TUBE	Trimline Two-Piece Reuse	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
2-BVD	BP PORT FITTING,2TUBE,ECON BLB	FlexiPort Accessory	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
2-BVDC	BP PORT FITTING,2 TUBE,ECONCLMP	FlexiPort Accessory	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
2-BVL	BP PORT FITTING,2-TUBE,LG BULB	FlexiPort Accessory	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)

[illegible]

Material #	Material Description	Cuff Style	NOTES ON CHANGE
30503514SA	CUFF, SOFT, LG AD, 1-TUBE, TP	FlexiPort SOFT	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
30506-121A	CUFF, VINYL, SM AD, 2-TUBE, BV	FlexiPort VINYL	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
30506-122A	CUFF, VINYL ADULT, 2-TUBE, BV	FlexiPort VINYL	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
30506122LA	CUFF, VINYL AD LONG 2-TUBE, BV	FlexiPort VINYL	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
30506-124A	CUFF, VINYL LG AD, 2-TUBE, BV	FlexiPort VINYL	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
30506-130A	CUFF, VINYL, AD, 2-TUBE, GAUGE	FlexiPort VINYL	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
30506-510A	CUFF, VINYL INFANT, 1-TUBE, TP	FlexiPort VINYL	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
30506-511A	CUFF, VINYL CHILD, 1-TUBE, TP	FlexiPort VINYL	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
30506-512A	CUFF, VINYL, SM AD, 1-TUBE, TP	FlexiPort VINYL	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
30506-513A	CUFF, VINYL ADULT, 1-TUBE, TP	FlexiPort VINYL	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
30506-514A	CUFF, VINYL LG AD, 1-TUBE, TP	FlexiPort VINYL	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
5082-03	CUFF, INFANT W/1 TUBE BAG	Two-Piece	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
5082-07	CUFF, NEWBORN, W/1 TUBE BAG	Two-Piece	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
5082-08	CUFF, INFANT W/2 TUBE BAG	Two-Piece	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
5082-21	INF SYST, CHILD CUFF	Two-Piece	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
5082-21H	INF SYST, CHILD, HOSPITAL PACK	Two-Piece	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
5082-22	INF SYST, ADULT CUFF	Two-Piece	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
5082-22H	INF SYST, ADULT, HOSPITAL PACK	Two-Piece	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
5082-23	INF SYST, LARGE ADULT CUFF	Two-Piece	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
5082-23H	INF SYST, LG ADULT HOSPITAL PK	Two-Piece	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
5082-24	CUFF, CHILD W/2 TUBE BAG	Two-Piece	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
5082-25	CUFF, ADULT W/2 TUBE BAG	Two-Piece	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)

Material #	Material Description	Cuff Style	NOTES ON CHANGE
5082-26	CUFF, LARGE ADULT W/2 TUBE BAG	Two-Piece	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
5082-42	CUFF, CHILD W/1 TUBE BAG	Two-Piece	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
5082-43	CUFF, ADULT W/1 TUBE BAG	Two-Piece	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
5082-44	CUFF, LARGE ADULT W/1 TUBE BAG	Two-Piece	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
5082-45	INF SYST, THIGH CUFF	Two-Piece	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
5082-45H	INF SYST, THIGH, HOSPITAL PACK	Two-Piece	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
5082-62	CUFF, CHILD PRINT, 1 TUBE BAG	Two-Piece	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
5082-77	CUFF, THIGH W/1 TUBE BAG	Two-Piece	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
5082-78	CUFF, THIGH W/2 TUBE BAG	Two-Piece	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
5089-01	INF BAG, ADULT 2-TUBE BAG	Two-Piece	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
5089-01H	INF BAG, ADLT 2-TUBES,12 / BOX	Two-Piece	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
5089-02	INF BAG, CHILD 2-TUBE BAG	Two-Piece	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
5089-02H	INF BAG, CHILD 2-TUBE, 12 / BOX	Two-Piece	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
5089-03	INF BAG, ADULT 1-TUBE	Two-Piece	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
5089-04	INF BAG, CHILD 1-TUBE	Two-Piece	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
5089-11	FEMALE TRI PURPOSE TWIST FIT	Two-Piece	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
5089-13	COIL TUBING, 8 FT LENGTH	Two-Piece	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
5089-13H	COIL TUBING, 8 FT 12/BX	Two-Piece	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
5089-18	INF BAG, THIGH 1-TUBE	Two-Piece	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
5089-19	INF BAG, THIGH 2-TUBE	Two-Piece	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
5089-19H	INF BAG, THIGH 2-TUBE 6 / BOX	Two-Piece	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
5089-21	INF BAG, LARGE ADULT-1 TUBE	Two-Piece	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)

[illegible]

[illegible]

[illegible]

[illegible]

[illegible]

[illegible]

Material #	Material Description	Cuff Style	NOTES ON CHANGE
VINYL-11L-2BVC	CUFF, VINYL AD LONG 2-TUBE, BVC	FlexiPort VINYL	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
VINYL-12-2BV	CUFF, VINYL LG AD, 2-TUBE, BV	FlexiPort VINYL	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
VINYL-12-2BVC	CUFF, VINYL LG AD, 2-TUBE, BVC	FlexiPort VINYL	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
VINYL-12L-2BV	CUFF VINYL LG AD LONG 2-TUBE BV	FlexiPort VINYL	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
VINYL-12L-2BVC	CUFF VINYL LG AD LONG 2-TUBE BVC	FlexiPort VINYL	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
VINYL-13-2BV	CUFF, VINYL THIGH, 2-TUBE, BV	FlexiPort VINYL	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)
VINYL-13-2BVC	CUFF, VINYL THIGH, 2-TUBE, BVC	FlexiPort VINYL	Tri-Purpose Connector can no longer function as a luer-slip or luer-lock; it can only function as a threaded connector that mates with the female tri-purpose (ex. 5082-100)